


User Instruction for the Basic Armchair Commode With 20L Buckets – Fixed Height T69

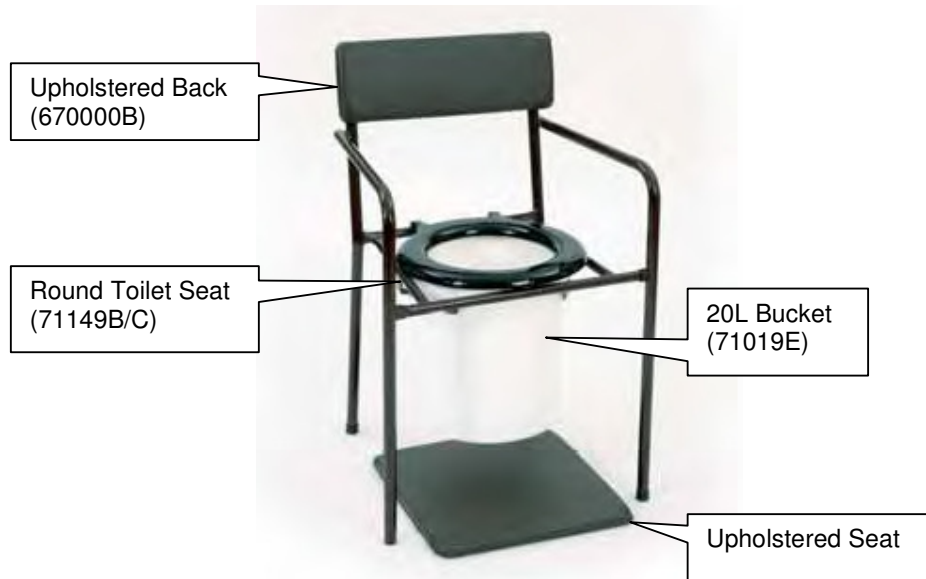

1. Introduction

This Basic Armchair Commode is designed with a raked back to enhance patient comfort during longer seating. It features heat sealed padded upholstery – this prevents water ingress. The height is adjustable and the feet are none-slip, none-marking for extra patient safety. The commode takes a capacity of 20 litres.

Product Dimensions	
Toilet seat height	520mm-680mm
Front width between arms	485mm
Back width between arms	419mm
Toilet seat to top of back	336mm
Toilet seat to top of arms	177mm
Toilet seat depth	356mm
Floor to arms height	692mm
Footprint	501mm x 540mm

2. Contents of the Box

- Basic Armchair Commode (T69)
- 20L Bucket (71019E)
- Upholstered Back (67000B)
- Upholstered Seat
- Round Toilet Seat (71149B/C)
- User Instructions (this leaflet)

3. Assembly Instructions

- Carefully unpack and examine the product to ensure it is complete.
- Before use check for any signs of damage e.g. cracked weld, flaking of paint, bent tube. If any signs occur please contact your supplier to exchange and refrain from using the product. **DO NOT ATTEMPT TO USE IF PRODUCT IS DAMAGED.**
- To remove the potty - lift the toilet seat and fit the plastic lid (to prevent spillage) before lifting the potty clear.
- To refit the potty - remove the plastic lid, replace the potty and then lower the toilet seat. If required, refit the vinyl seat.

4. Correct Use

User Instruction for the Basic Armchair Commode With 20L Buckets – Fixed Height T69

- Please empty the potty at regular intervals and be aware of the weight limit.
- Do not attempt to move the commode whilst occupied.
- Remove the vinyl seat when you want to use the chair as a commode.

5. Cleaning

- Use mild detergent with warm water to wipe down framework on at least a weekly basis more frequently if required.
- Pay particular attention to the hinges and the underside of the toilet seat.
- Do not use an abrasive cleaner as this will damage the frame and nullify your warranty.

Use of power/jet washing & industrial washing machines:

- If either of the above is used to clean equipment care must be taken since high-pressure water can be penetrative.
- Care should be taken that frames are not placed too close to high-pressure water jets.
- After cleaning ensure any residual water is removed from the framework as it could cause permanent rusting of the frame (leave frames open in a warm environment for several hours to ensure all internal moisture evaporates).
- We cannot offer advice on industrial washing machines with regards to specific machine settings, as these will be dependant on the manufacturer's instructions.

6. Maintenance Requirements

The product should be checked on a regular basis for damage or wear and tear – this includes flaking or cracking paint. Any undue signs of damage should be noted and reported to the supplier in the first instance. Should damage be evident, then the user should refrain from using the product until advice has been sought from the supplier.

7. Warning

- This product is intended for users of a maximum of 125kg/20stone in weight.
- We recommend a full risk assessment should be carried out before use.
- Under no circumstance should the commode be moved whilst occupied.
- It is advised that the frame is checked periodically for signs of wear and tear. Any problems please report to your supplier immediately.

8. Warranty

This product carries a 12 month warranty from the date of purchase against manufacturing defects. Misuse or damage is not covered. This warranty does not affect your statutory rights. As manufacturers we cannot specify expected product life as this is based on the level of usage.

9. Disposal Instructions

- Appropriate hygiene standards should be applied before disposal to reduce the risk of infection.
- Materials used are essentially non hazardous and can be disposed of appropriately.

10. Contact Details

We have attempted to cover most eventualities, however if you have any specific queries, please contact:

Cefndy Healthcare
Cefndy Road

Rhyl
Denbighshire
LL18 2HG

Tel: 0044 1745 343877 (Quality Assurance Dept)

E-mail: cefndy.sales@denbighshire.gov.uk